


# Fluor *i* Gel

**DNA gel**

**LED light**

**Fluorescent detection**

Fluor *i* Gel is a device that images DNA electrophoresis gel and analyzes. Fluor *i* Gel Blue is optimized for the reagents, wavelength of 450nm~490nm, that is invented alternative to EtBr. It consists of LED light and dark room chamber. Fluor *i* Gel can also be used in conjunction with UV transilluminator. Compact size and simple design has satisfied user's convenience. Gel can be observed through the window at the top and gel cutting can be done conveniently through the doors at each side. Our simple program enables users to obtain results easily and simply quantify the DNA band.

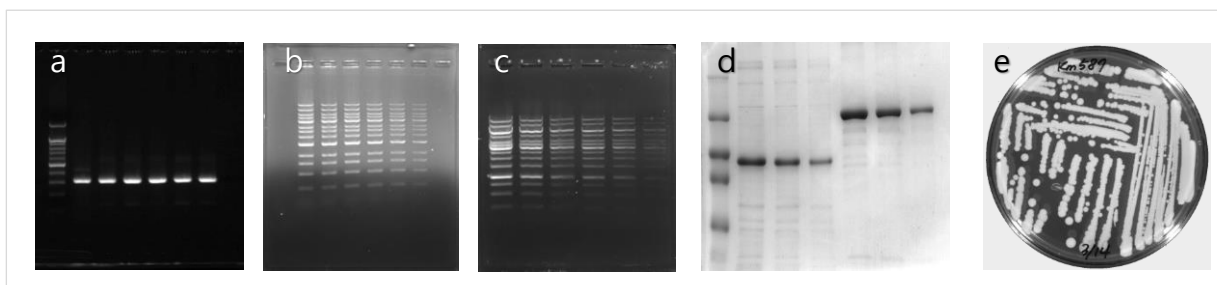

a. Plasmid DNA stained with loading type EtBr alternative reagent. b. DNA marker electrophoresed on agarose gel containing a gel mixture type of EtBr replacement reagent. c. DNA marker electrophoresed on agarose gel containing EtBr. Images obtained using a UV transilluminator and Fluor *i* Gel in combination. d. PAGE gel image. White plate, white light used. e. Petri dish image. White plate, white light used.

## Specifications

Camera	1/2" 1.3M 8bit CMOS, 1280 x 1024 pixels
Size (W x D x H)	260 x 260 x 400 mm
Interface connector	Standard USB 2.0
Field of View	160 x 130 mm
Light source	Blue LED
Software	Image capture, Set ROI: Manual or Automatic
Analysis	Subtract background, measuring of intensity

## NEOgreen

DNA staining reagent  
(EtBr alternative reagent)

NEOgreen is developed with the purpose of replacing EtBr. Same as the existing method, to make 100ml of Gel, 5~10ul of NEOgreen is needed. Therefore, the experimental method doesn't need to be changed. Blue light (470nm) and UV transilluminator can be used with it to observe the DNA gel



## NeoScience Co., Ltd.

#111 Toyou 2cha 131 Munwhawon-ro Yuseong-gu, Daejeon Republic of Korea 34184  
(+82) 70-7430-6829 [www.invivoimage.com](http://www.invivoimage.com) [neoscience2011@gmail.com](mailto:neoscience2011@gmail.com)